

# SKYMEN Spare Parts Benchtop Ultrasonic Cleaner 6.5L 30min Timer

#### **Basic Information**

Place of Origin: ShenzhenBrand Name: Skymen

Certification: CE/RoHS/ISO9001

Model Number: JP-031Minimum Order Quantity: 1

Price: usd 140-175/pcs
Packaging Details: Wooden Case
Delivery Time: 5-8 days

Payment Terms: T/T, Western UnionSupply Ability: 3000 PCS per month



## **Product Specification**

• Tank Size: 300\*150\*150mm

Tank Capacity: 6.5L
Ultrasonic Power: 180W
Timer: 0-30min
Transducer: 3 Pcs
Frequency: 40KHz

Heater: 20-80 Degrees AdjustableN.W/G.W: N.W.5.4KG G.W.6.2KG

• Highlight: Spare Parts Benchtop Ultrasonic Cleaner,

30min Timer Benchtop Ultrasonic Cleaner,

6.5L Lab Ultrasonic Cleaner



## **Product Description**

Fully Clean Dirt Strong Power Quick Delivery Pump Benchtop Ultrasonic Cleaning Machine Ultrasonic Cleaner Bath



		Drainage	no
5		Cover	include
		Spare part	Basket, latex gloves, warranty card, certificate of quality, manual
	Pa cki ng Fe atu res	Tank Size	300*150*150mm( L x W x H )
		Unit Size	375*200*275mm ( L x W x H )
		Package Size	410*260*350mm ( L x W x H )
		N.W/G.W	
		Warranty	1 year
		2. Power adjustable by mechanical panel 3. Stainless steel washing basket and Lid, 4. ultrasonic cleaning solutions penetrate virtually any soil type including: oils, residues, solvents, pastes and chemicals and other contaminants are safely lifted and dispersed, leaving a surface appearance that is clean and free of particulate; even in crevasses and internal bores.	
6	Ap pli cat ion	Ultrasonic clean, emulsify, mixing, extraction, disperse	
7	Ce rtifi cat es	CE, ROHS, FCC, ISO9001:2008	

#### **Special Features:**

Digital control Timer and Heating Temperature control function

Isolated-Circuit Design, Keep Water Away for Extra Safety

Powerful Transducer for Best Cleaning Result


Top quality steel made for body, cover and base

Basket and Tray Anti-noisy and Shock absorption Handle Drain Sophisticated made and customized design Professional testing, quality assurance

### Application:

- -Scientific Labs:Eyeglass Frames;Lab Glassware;Optical&Contact Lenses;Pipettesscientific Instruments;Test tubes
- -Industrial Manufacturing: Assemblies; Gears; Metal & Plastic Parts; Precision Bearings
- -Electronic manufacturing:Capacitors;Ceramic Substrates;PC Boards;Packaging Components;Quartz Crystals
- -Medical and Dental Lab:Burs;Blood Oxygenators;Dental Instruments;Syringe Parts;Surgical Instruments
- $\hbox{-Jewelry Manufacturing:} Chains; Charms; Coins; Intricate Settings; Precious Metals \& Gemstones; Watches$

If your company has other needs, you can choose the size that suits you according to the chart below

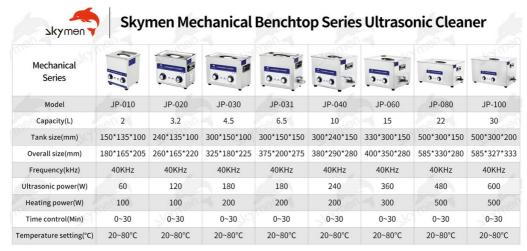

#### How ultrasonic cleaner works:

Electrical energy is converted into mechanical vibrations thanks to special elaborated transducers. These vibrations are transferred into the fluid through the walls of stainless steel. This creates tiny microscopic vacuum bubbles, which cause consecutive implosions (cavitation). The high-energy jets that are thus created can efficiently remove dirt from surfaces an materials placed in the cleaning bath.

An ultrasonic cleaner in your kitchen will save time and reduce effort while providing a safe solution for your kitcheware. With the help of oil remove detergent ,the ultrasonic cleaner could fast and efficient take the oil ,polish the kitcheware.

when using an ultrasonic cleaner on glassware, plasticware, and your lab instruments. Whether you are cleaning oven trays, pans. studies have shown that it can be up to 30 times faster while delivering thoroughly clean results.



#### **Skymen Technology Corporation Limited**



86-755-27094405





skymenultrasonic.com

Floor 1st & 2nd, Building 3, Tanggang Taifeng Industrial Park, Dawangshan Community, Shajing Street, Bao'an District, Shenzhen